

# FERNO Compact 2 Track

Carry chair and track chair combined in one –  
the ultimate in versatility and easy handling



FOR MORE DETAILED  
INFORMATION SEE OVER

## Specifications

Height	Depth	Width	Weight	Load Limit
810 mm (folded)	295 mm (folded)	450 mm	15 kg (chair+track)	200 kg/31st/434lbs
1062 mm (deployed)	950 mm (track deployed)		11.4 kg (chair)	
1445 mm (rear handle fully extended)			3.6 kg (track)	
485 mm (seat height)			Ferno storage bracket 7.5 kg	
Track module length	Track module width	Track module depth		
730 mm	340 mm	120 mm		

**FERNO** – is a world leading innovator and manufacturer of patient transportation solutions.

Our product range is born of a closer understanding of our customers' needs and the practical demands of their work. Ferno's industry-leading **Compact 2 Track Chair** typifies our passion for top quality, highly functional design... patient transfer products that don't just do the job, they do it **BETTER**.

call +44 (0) 1274 851 999  
email [sales.uk@ferno.com](mailto:sales.uk@ferno.com)  
www [ferno.com](http://ferno.com)



# FERNO Compact 2 Track




**A lightweight, manoeuvrable carry chair... with the capability of converting into a track chair for downstairs transfers**

Having all the advantages of the familiar and commonly used **Compact 2 carry chair**, lightweight, easy to handle and robust, the **Compact 2 Track** also has the capability of being able to attach a track to the chair as and when needed, thereby allowing patients to be safely and easily transferred downstairs, eradicating the need to lift the chair occupant and thereby minimising operator injury risk.

An adjustable-height rear handle also ensures that the chair accommodates operators of all statures, eliminating stooping when moving patients on vertical ground, up and down vehicle ramps etc, making the transfer easier, safer and more stable. This adjustability and improved operator stance also gives more stability to the chair when traveling over rough ground.

**Combining the features of a carry chair and track chair in one product reduces manual handling risks, reduces the weight of equipment operators have to carry, decreases the equipment weight in the vehicle and frees-up available space, which all ultimately lead to cost savings.**

- **Outstanding versatility** – converts to a track chair quickly and easily, as and when needed
- **Enhanced operator use** – rear handle adjusts for operators of different stature
- **Quick and easy to use** – Simple chair deployment, rear handle height adjustment, track attachment and removable headrest
- **Lightweight and durable** – stainless steel construction, designed to minimise weight
- **Infection Control** – seat cover and headpad made from materials containing anti-bacterial additives
- **Stays stable, even on uneven ground** – adjustable rear handle, large wheels and axle bar ensure excellent balance
- **Minimised injury risk** – reduced weight of equipment to carry, improved operator body posture, elimination of lifting of patients downstairs and in/out of ambulance vehicles
- **Easy storage** – in-vehicle 2 piece storage bracket secures both the chair and track in one place ready to be deployed or, store as a standard chair keeping the track separately in a carry bag

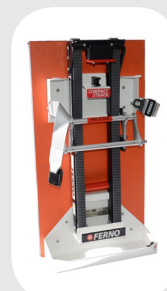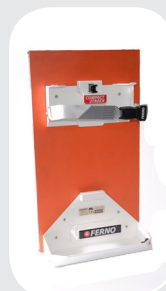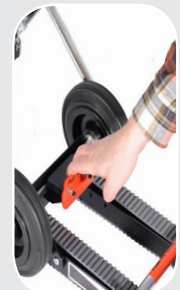

2-Piece In-Vehicle Storage Bracket for chair & track

**Part Nos:**

071481450  
009103800

Compact 2 Track Chair  
Compact 2 Track In-Vehicle Storage Bracket for chair & track

FERNO - Compact 2 Track FPC010 03.25

**call** +44 (0) 1274 851 999  
**email** sales.uk@ferno.com  
**www** ferno.com



Ferno (UK) Limited, Stubbs Beck Lane, Cleckheaton, West Yorkshire, BD19 4TZ, England