

POWERXTM MK2

Frontline & PTS Ambulance Stretcher



The POWERX MkII enhances the form and functionality of X-frame Ambulance stretchers, ideal for use with today's tail-lifts and ramps.

Feedback from end-users highlighted the need for a more robust stretcher tailored for the pre-hospital market. Its ergonomic design and enhanced features improves operation, handling and manoeuvrability, making it safer for operators while enhancing the patient comfort experience.

POWERX MkII is one of the most configurable and feature-rich Ambulance Stretchers available, placing the latest technology at your fingertips.



Push-Poles with ball handles



Latching Backrest Defib Hook



Simple to operate controls



Easy to apply brakes and directional wheel-lock



Optional CD/ZA 02 cylinder holder

POWERX[®] MK2

Features


- **Power that is dependable and reliable** – proven lithium-ion battery technology, quick to recharge and long-lasting during daily operation and throughout the stretcher's lifetime
- **Range of Features and Options** – usable over a range of differing needs and purposes, both Frontline and PTS with future-proofing add-ons
- **Quick configurability** – integrated side-rail interface on each side for accessories to be easily attached, making the stretcher more immediate to differing clinical needs
- **Easily re-configured** – A&E to Specialist Care and Bariatric, all from one single stretcher design
- **Integrated infection control** - clean lines, wipe-clean pressure-reducing AB-PR mattress, biosafe patient restraints and anti-bacterial paint surfaces all contribute to a high level of safety
- **Simple controls** – up and down at the touch of a button. Optional side lights make it highly visible in poor conditions
- **CEN Tested Safety** – dynamically tested in six directions and certified to BS EN 1789

POWERX comes with the most advanced choices of features and configurable options ever developed in a patient transport system

- Smooth powered raise/lower movement, utilising latest electro-mechanical actuator
- Robust design actuator for exceptional stability in the raised position
- Unassisted lift of 350 kgs (55 st) SWL up to the maximum raised trolley height of 880 mm, 400 kgs (63 st) in lowered position
- Sturdy head-end & foot-end Push-poles
- MCF winch-link connection for assisted loading up vehicle ramps
- Anti-Bacterial Pressure-Reducing mattress
- Wipe-clean Biosafe patient harness
- Adjustable winged head-rest (AB-PR)
- Wheel brakes on all 4 corners
- Directional wheel locks at head-end
- SX Surface Extending cot-sides
- Side-trolley lighting option
- Latching backrest monitor hook option
- CD/ZA O2 cylinder holder option
- Supplied with one Milwaukee M28 5Ah battery as standard



Available Accessories

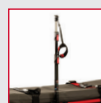

3-Stage IV Pole quickly place an IV anywhere around the patient
LAIV3SCL-UK




SX Surface Extender provides added comfort and security for patients of any size
FWESX



PacRac+ Equipment Mount with iNTRAXX profile supports monitoring devices, opening both sides
FWEPR



Storage Net keep items to hand when transferring patient to and from the vehicle
082-2453

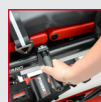

Side flip-out handles, aids manoeuvrability in a range of environments
FWESLH




Split-Scoop bariatric kit giving full length patient surface width extension for bariatric patients
009008250

Split-Scoop foam mattress sidepads
0450-1023 left, 0450-1024 right



Spare Milwaukee M28 5Ah battery
75-0710-001
240V AC and or 12V DC battery charger
75-0711-002

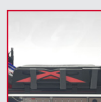

SX Cot-side Extender Wings, giving more surface area for larger patients

Specifications

LENGTH	Shortened / Extended	1650 / 1930 mm
HEIGHT	Low to Raised	480 / 880 mm
WIDTH	Normal	570 mm
	SX Cot-sides horizontal	920 mm
	Split Scoop kit attached	945 mm
Load Limit (SWL) Unassisted lifting	Raised	350 kg / 55 st
	Lowered	400 kg / 63 st
POWER	28 Volt DC Battery 5Ah	
WEIGHT		From 79 kg*
FASTENER	Two-Part Lock	13 kg

* Weight dependant on final configured specification. POWERX illustrated shown with various optional configurations and accessories. For full available range see Configurable Options and Accessories list.