

FERNO HARRIER LT-LBS

Bariatric Ambulance Trolley

Avoid operator injury with this dedicated bariatric powered trolley that lifts up to 350 kg unassisted and with extending cotsides that give increased patient surface as and when required



FOR MORE DETAILED
INFORMATION SEE OVER

Specifications

Height	Height	Length	Length	Width	Width LBS	Weight	Lift Load	Load
Raised	lowered		Folded head end		Extended		Limit	Limit
							Raised position	Lowest position
925 mm	480 mm	1930 mm	1610 mm	570 mm	760/955 mm	From 69 kg*	350 kg	400 kg
						*Dependant on configuration	55st	62st
							770lbs	868lbs

FERNO – is a world leading innovator and manufacturer of patient transportation solutions.

Our product range is born of a closer understanding of our customers' needs and the practical demands of their work. **The Harrier LT-LBS Trolley** typifies our passion for top quality, highly functional design... patient transfer products that don't just do the job, they do it **BETTER**.

call +44 (0) 1274 851 999
email sales@ferno.co.uk
www ferno.co.uk



FERNO HARRIER LT-LBS

Bariatric Ambulance Trolley

Exceptional lifting power – with up to 77 lifts from a single battery charge




Ambulance trolleys need to be easily manoeuvrable and quickly secured within the vehicle – these considerations are an absolute must. But to perform the task to maximum efficiency you also need a trolley which can lift and carry even heavier patients without risk of injury to the operator – and without needing frequent battery recharges.

Ferno's **HARRIER LT-LBS** combines exceptional lifting power with a high capacity battery. It is capable of lifting up to 350kg unassisted, and will perform an amazing 77 lifts of 100kg with just a single battery charge.

The **HARRIER LT-LBS** is configurable to differing customer needs with a variety of options and accessories available, designed to improve the ease of use of the trolley.

- **Reduction in manual handling injuries** – lightweight, easy to operate design, incorporating ergonomic push/pull handle options
- **Smooth powered lifting action** – lift up to 350 kg, with a soft start/stop vertical movement
- **Easy to manoeuvre** – lightweight trolley, using optimum performance low resistance wheels, with various operator push/pull handles and push-pole options
- **Simplified high-performance battery management** – Removable and rechargeable 28 V battery, with up to 77 lifts from one charge, visible battery level indicator, and just 1 hour recharge time (2 batteries & 240V charger supplied)

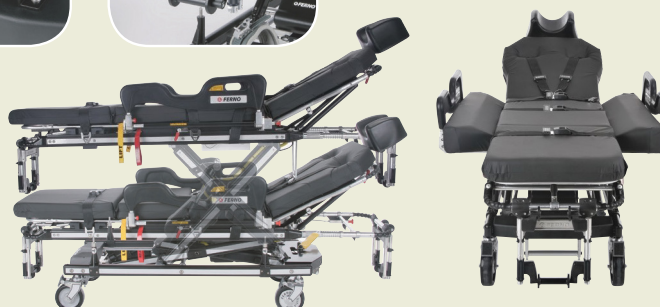

- **Excellent hygiene and infection control** – IPC compliant, incorporating anti-microbial paint finish, use of Silver Biosafe harness restraints and easy to clean pressure-relief mattress
- **Safe at night** – Reflective X-frame for increased visibility in low levels of light

Trolley configurations:

Harrier LT-LBS, standard-length extending cotsides (730 mm) with padded sidepads (900 mm)
Harrier LT-LBS, long extending cotsides (895 mm) with inflatable sidepads (1400 mm)

*Both configurations, cotsides extend to 2 width positions, 760 mm and 955 mm width.

Full range of configurable options and accessories available:

Directional wheel lock
Anti-static wheel
Detachable vertical push-poles at foot-end
Folding push-poles at foot-end (not available on LT-LBS with long cotsides)
Push/pull handle at foot-end
Removable adjustable winged headrest pad
IV Pole (detachable)
Black webbing patient restraints instead of Biosafe
Winch link attachment
CD cylinder holder (on foot-end)
PacRac equipment holder
12V battery charger
Patient shield
Trolley dust cover, grey AB material

Available with various locks: Standard 2-part lock, NMI Track lock, LHS-operated NMI Track Lock, Unwin Track Lock, LHS-operated Unwin Track Lock, SlideLock 150, SlideLock 255

Bundle part number dependant on configuration and lock required

FERNO - Harrier LT FPC005 08.21

call +44 (0) 1274 851 999
email sales@ferno.co.uk
www ferno.co.uk



Ferno (UK) Limited, Stubs Beck Lane, Cleckheaton, West Yorkshire, BD19 4TZ, England