


# **FERNO** Neil Robertson Stretcher

Industry standard stretcher ensuring safety during confined space rescue and vertical or horizontal lifting

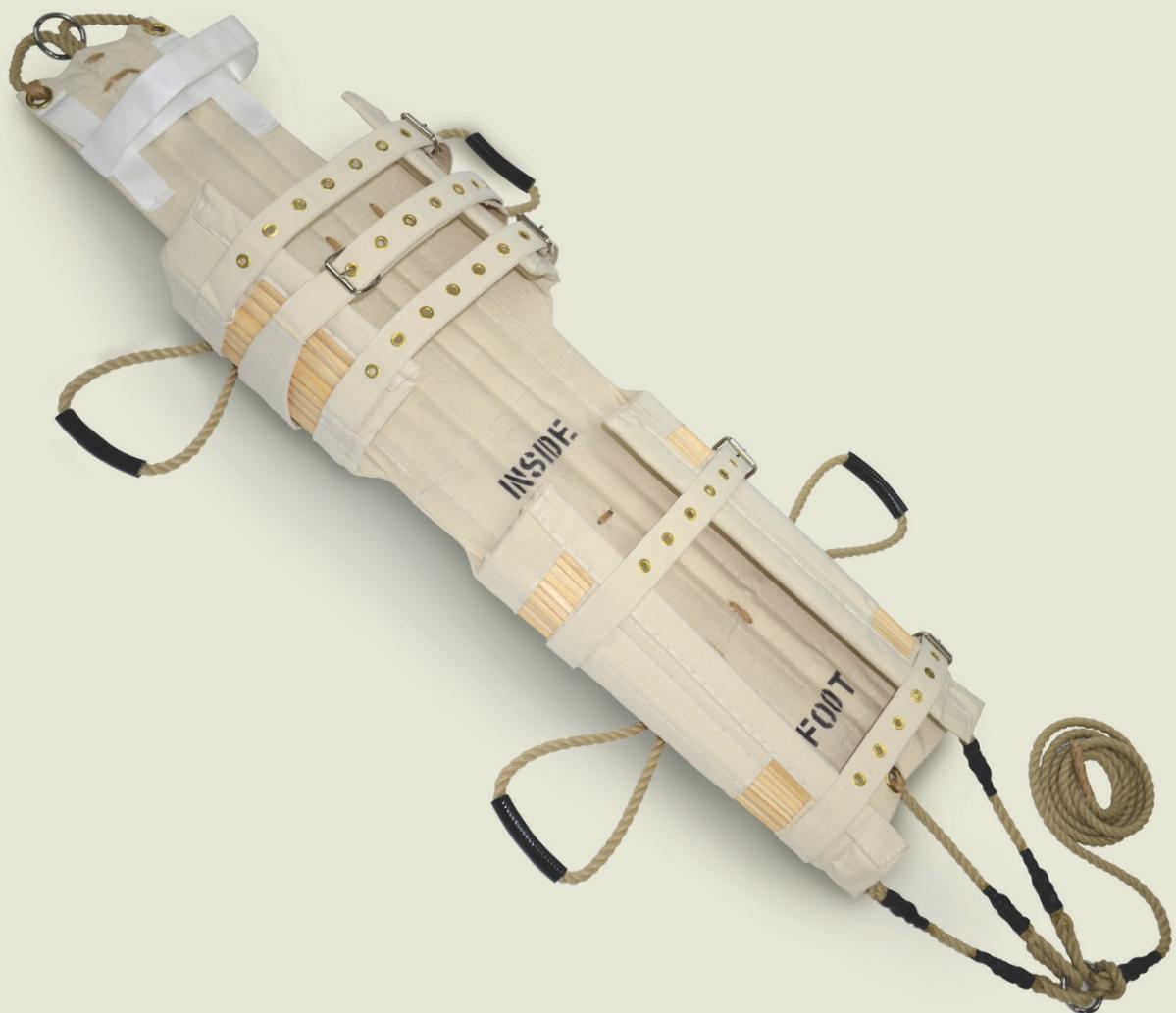

**FOR MORE DETAILED  
INFORMATION SEE OVER**

## Specifications

**FERNO** – is a world leading innovator and manufacturer of patient transportation and rescue solutions.

Our product range is born of a closer understanding of our customers' needs and the practical demands of their work. Ferno's industry-leading **Neil Robertson Stretcher** typifies our passion for top quality, highly functional design... patient transfer products that don't just do the job, they do it **BETTER**.

### Neil Robertson Stretcher

Length	Width (open/folded)	Depth (open/folded)	Weight	Load Limit
1500 mm	1000 mm/310 mm	40 mm/180 mm	8 kg	136 kg/21st/300lbs

**call** +44 (0) 1274 851 999  
**www** ferno.co.uk

**email** sales@ferno.co.uk



# **FERNO** Neil Robertson Stretcher

## Industry standard stretcher to ensure stability before and during lifting

For many years the Neil Robertson's unique 'spliced-slat' design has been the stretcher of choice in difficult rescue situations, where a casualty needs to be manoeuvred, lowered or lifted to safety. It is widely used by the Royal Navy and Merchant Navy, in Mining and Construction Industries, and by a number of Emergency Rescue Services.

- **Extremely tough** – made from spliced, ram-in slats with a rot-proof cotton exterior
- **Strong and supportive** – retaining straps, under-arm flaps and built-in head restraints offer stability during vertical lifting
- **Safety that has stood the test of time** – an acknowledged standard used in industry
- **Enhanced confidence** – durable manilla support ropes, carry handles, hoisting rope and guide rope ensure security when moving

## Accessories available:

- Carrying bag for convenient storage when not in use
- Pillow for head or neck support

## Storage/carry bag



## Head/neck support pillow



## Part Numbers:


010779910 Neil Robertson Stretcher, 1.5m

010779920 Neil Robertson Stretcher, 1.5m with neck/head pillow

## Accessories:

031347500 Carrying Bag for Neil Robertson 1.5m

010779940 Neck/head pillow for Neil Robertson

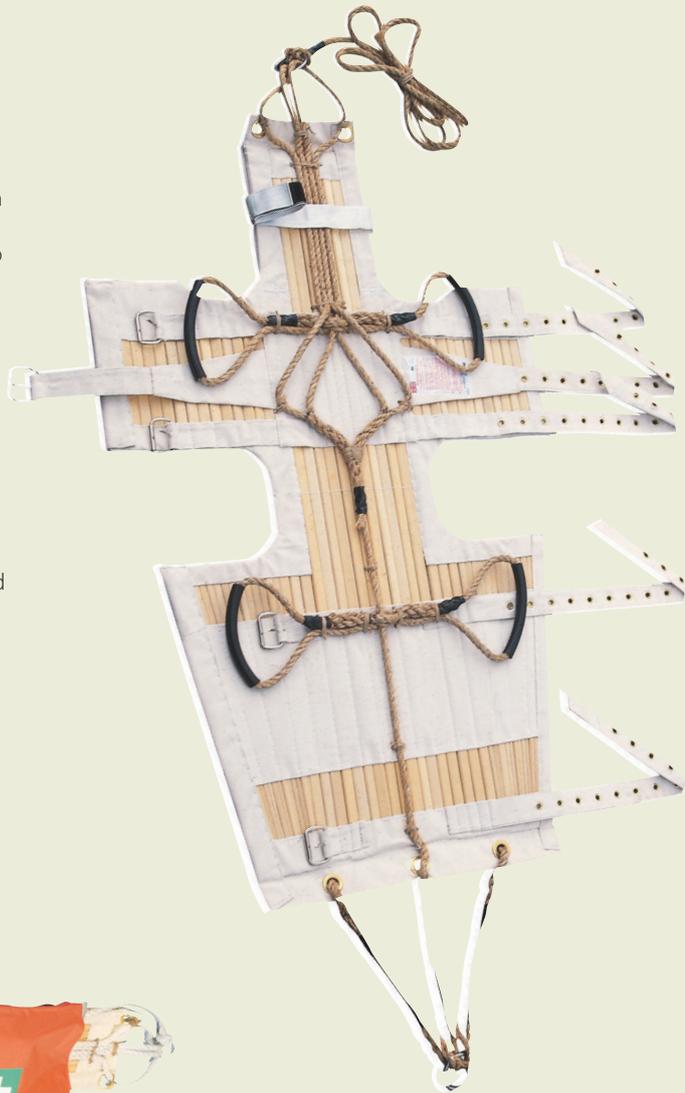

Vertical lift using integral lifting O-ring and foot-end guide rope

FERNO - Neil Robertson Stretcher FPC006 0320

**call** +44 (0) 1274 851 999 **email** sales@ferno.co.uk  
**www** ferno.co.uk

Ferno (UK) Limited, Ferno House, Stubbs Beck Lane, Cleckheaton, West Yorkshire, BD19 4TZ England.

